

Real patients. *Real results.*



Before and after 1 treatment



Courtesy of Zena Gabriel, MD



Before and after 1 treatment



Courtesy of Fine Fittings Med Spa



Before and after 1 treatment



Courtesy of Melanie Palm, MD



Before and after 1 treatment



Courtesy of Fine Fittings Med Spa

## Follow us

 [cynosurelasers](#)

 [Cynosure, LLC](#)

 [lutronic\\_us](#)

 [Lutronic](#)

[www.cynosurelutronic.com](http://www.cynosurelutronic.com)

## Frequently asked questions

### Is it painful?

Everyone's experience is different. Your provider will choose the best way to keep you comfortable during treatment. Usually, a numbing cream is applied to help with treatment comfort.

### What is the downtime like?

Downtime depends on the treatment:

**Precision Peel:** Mild redness for 6 to 48 hours.

**Deeper resurfacing:** Redness can last a week or longer. Makeup should be avoided for two days. Most people can return to regular activities after about 5 days. Your provider will explain what to expect based on your specific treatment.

### What results can I expect?

The results after a single medium to aggressive treatment can be dramatic, and skin improvement benefits can extend well beyond 6 months. For lighter treatments, the skin improvement effect may be more subtle. Your provider can help with the proper treatment plan to meet your goals.

# eCO2<sup>TM</sup> 3D

## Reveal the Real You



## CYNOSURE® LUTRONIC®

Intended for distribution in US and Canada.

The eCO2 3D aesthetic workstation is indicated for the following: Skin resurfacing, wrinkles and fine lines, furrows, scars, acne scars, surgical scars, burn scars, textural irregularities, dyschromia, pigmented lesions, lentigines, actinic & solar elastosis, actinic keratosis, precision ablation and/or vaporization of soft tissue in dermatology and plastic surgery, coagulation & resurfacing of soft tissue, keloids, xanthelasma palpebrarum, syringoma, laser derm ablation, debridement, and vaporization, ablation, vaporization, excision, incision, and coagulation of soft tissue, and rhinophyma. Cynosure, LLC and Lutronic Aesthetic, Inc. do not diagnose or treat medical conditions, prescribe treatment protocols, provide medical advice, or otherwise practice medicine. Patient should consult with a qualified healthcare professional for diagnosis and proper treatment protocols/parameters to determine if treatment with any product is right for them. Qualified healthcare professionals are responsible for diagnosis and treatment protocols using products.

© 2025 CYNOSURE LUTRONIC, INC. All Rights Reserved. eCO2 3D, Precision Peel, and Controlled Chaos Technology are trademarks and Lutronic is a registered trademark of Lutronic Corporation. Cynosure is a registered trademark of Cynosure, LLC. PRD-2059

921-7077-003 R02 US/CAN-EN

# Your skin, the way you want it

The eCO2 3D workstation is the most versatile CO2 laser in aesthetics, designed for gentle touchups, total transformations, and everything in between. You can choose gentle, low-downtime sessions with **Precision Peel** or a single treatment with more dramatic results. Your provider can help you create a treatment plan that's right for you.


## Common results:

- **Softens** fine lines and wrinkles
- **Improves** skin texture and the appearance of scars
- **Targets** pigmentation and vascular conditions
- **Reveals** bright and firm skin



## How it works

eCO2 3D treatments help improve skin by using your body's natural healing process to create new tissue and collagen, making your skin look younger and healthier.

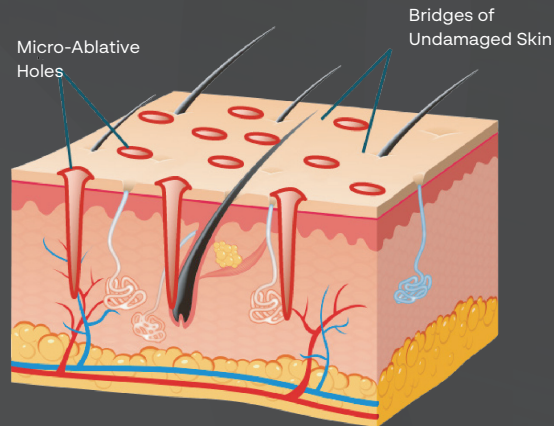


## Treated skin histology

1-2mm deep



### Immediately after treatment

Fractional laser energy creates micro-ablated columns to stimulate cellular regrowth. The surrounding healthy tissue promotes rapid healing.

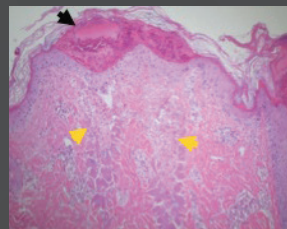

### Two days post-treatment

Skin is healing on the surface and creating new dermal collagen. This will continue for up to six months. The improved skin quality will make for a more youthful appearance.

## PRECISION PEEL by eCO2 3D



## Resurfacing redefined

**Precision Peel** is eCO2 3D's signature low-downtime resurfacing protocol, designed for patients seeking visible improvement in texture, tone, and radiance without the recovery associated with traditional CO2 treatments. Powered by the system's three unique tips, **Precision Peel** treatments can be performed quickly and comfortably, making them ideal for new patients, maintenance regimens, or event prep.